



OASYS / OASYS Hinge Research & Development Lab
660 Commerce Drive, Ste. C * Roseville, CA 95678
(916) 865-4528 www.dreamsystemsdentallab.com



510K—#K971794

OASYS Simplicity

Simple Mandibular Advancement
Unattached, Two Piece Appliance
Labial Shield for Advancement
Upper 2mm Clear Cushion
(Waxed Facial #6-11 for Reduced Labial Pressure)



EMA ThermoFlex Hard

Made With or Without Clasps
Resin Compensates of Impression
Defects, Missing Teeth
Flat Plane Opposing Splint Design



OASYS Standard Design

Adjustable, Intra-Oral Nasal Dilators
30%+ Patients Could Benefit
Improved Breathing Throat & Nose

EMA VacuForm Splints

Standard with 2mm Splints
Available in 2.5 MM Splints for
Bruxers, or in Hard / Soft Splints
Requires Very Accurate
Impressions or Models



OASYS Supreme

Improved Breathing Throat &
Nose, with Tongue Repositioning
Sliding Advancement with Locks
MM Scale...8-15mm Advancement



EMA First Step

90 Day Trial Device
0.75 MM Vacuform Splints
Three Yellow Strap Sizes 21,19,17mm
Requires Very Accurate
Impressions or Models



The O2 OASYS Oral / Nasal Airway System Uniqueness with Many Applications

- Simple Mandibular Advancement
- Advancement to 15mm
- Nasal Obstruction / Resistance
- Severe OSA
- Upper Denture
- Minimal Dentition Upper / Lower—Easily Modified to Function
- In Combination with Invisalign

510K—#K030440



Advancement Strap Pack

Straps provided for 5mm of
Advancement from Initial Setting
5 Pairs Blue, 2 Yellow
Bites Sent Advance of End to End
will Receive 1-2 Pairs of Straps to
Reduce Advancement, if Needed.

Available Options

Nasal Dilators
Tongue Lifter Buttons

Straps

Sizes

21mm

20mm

19mm

18mm

17mm

16mm

15mm

14mm

13mm

White—Very Flexible

(Provided Only Upon

Request)

Yellow—Flexible

Blue—Standard

Strap Less Flexible

Clear—Rigid

Provided Only Upon

Request for Less

Movement Side to Side





OASYS Research & Development Lab


Premium & Medicare OASYS Herbst Designs

Premium Welded

Fully-Welded Frame, OASYS Herbst Hinges with Clasps, Ball Hooks & Elastics


Medicare

OASYS Herbst Hinges with Welded Frame, Clasps, Ball Hooks & Elastics

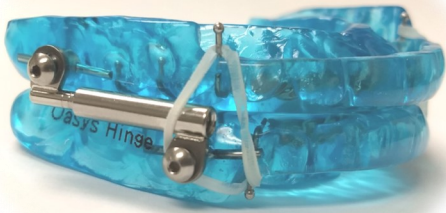

10 MM Advancement

MM Scale Shows Advancement
Set Open 1-2mm for Retraction
Shorter Length for Better Placement

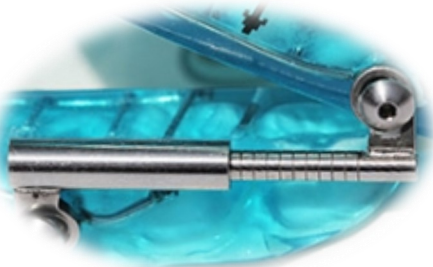


Intra-Oral Activation

Easy Advancement

One Revolution = 1/2 MM

Two Revolutions = 1 MM

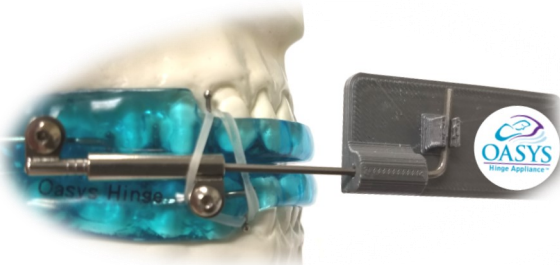

Available Options

Nasal Dilators
Tongue Lifter Buttons
Comfort Caps Available for
Side Screw Irritation
DentiTrac Compliance Monitor



510K—#K083209
Medicare Compliant
#M001 OASYS Hinge Appliance
+Nasal Dilators—#M002
+Tongue Lifters—#M003



Premium & Medicare Herbst Designs

Premium Welded

Fully-Welded Frame, Lingual & Buccal, with Herbst Hinges, Clasps, Ball Hooks & Elastics

Medicare

Welded Frame with Herbst Hinges
Clasps, Ball Hooks & Elastics



5 MM Advancement

Directional Arrows

Right Screw Up—Left Screw Down
Four Revolutions = 1MM—16 Clicks



Additional Advancement

2 MM Crimpable Shims Available
Longer Arms for Additional 5 MM
Comfort Caps Available for
Side Screw Irritation
DentiTrac Compliance Monitor



510K—#K083209
Medicare Compliant—#MC1
+Nasal Dilators—#M002
+Tongue Lifters—#M003

Invisalign Retainers



Upper Invisalign Tray
Used in place of Clear Aligner



**The Lower OASYS
Splint is Designed to
Function with Lower
Tray Transitions**

**“Reducing Air
Pressure Needed
Up to 40%!”**


**Improve Results -
Oral Appliance
with CPAP Therapy**



**“The CPAP Pro Bracket is
easily attached for Home
Use or Removed for Travel”**



Minimal Dentition

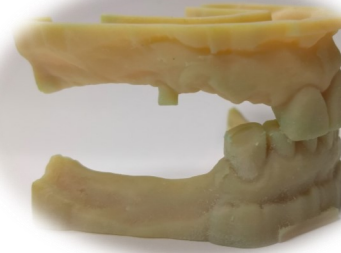


Maxillary Denture



OASYS
TMJ SPLINT
TONGUE TRAINER
SPEECH THERAPY



**Digital Scan /
Model Service**

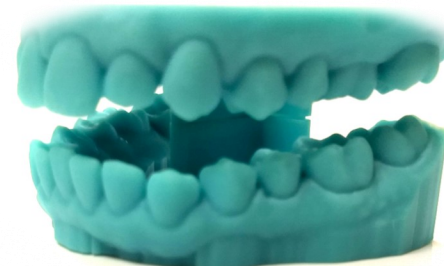

AM Aligner
Morning Repositioner
Anterior or Full Arch Thermo Wafer
for Fabrication on Patient or Preformed

