



## Additional Sleep Appliance Options



OASYS / OASYS Hinge Research & Development Lab  
660 Commerce Drive, Ste. C \* Roseville, CA 95678  
(916) 865-4528 [www.dreamsystemsdentallab.com](http://www.dreamsystemsdentallab.com)

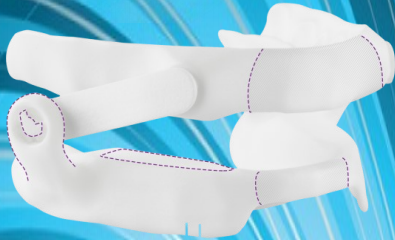

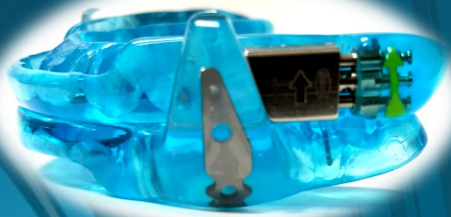

**OASYS Simplicity**

Simple Mandibular Advancement  
 Unattached, Two Piece Appliance  
 Labial Shield for Advancement  
 Upper 2mm Clear Cushion  
 (Waxed Facial #6-11 for Reduced Labial Pressure)



**OASYS Standard Design**

Adjustable, Intra-Oral Nasal Dilators  
 ≥30% of Patients could Benefit from  
 Improved Breathing in Throat & Nose



**OASYS Supreme**

Improved Breathing Throat &  
 Nose, with Tongue Repositioning  
 Sliding Advancement with Locks  
 MM Scale...8-15mm Advancement



**510K—#K030440**



**The O2 OASYS  
 Oral / Nasal Airway System**  
**Uniqueness with Many Applications**

- Simple Mandibular Advancement
- Advancement to 15mm
- Nasal Obstruction / Resistance
- Severe OSA
- Upper Denture
- Minimal Dentition Upper / Lower—  
Easily Modified to Function
- In Combination with Invisalign




**MyoSplint—Pivotal  
 Splint with Tongue Guides**



**Pivotal Splint  
 on First Molars**




**Pedo Speech  
 Trainer**

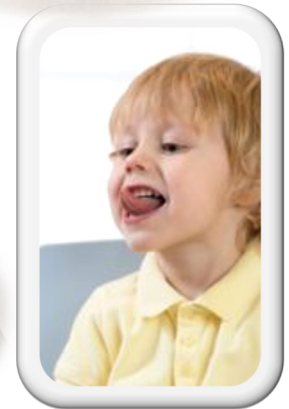

**Digital Scan / Model Service**

- ♦ iTero
- ♦ Trios,
- ♦ Carestream
- ♦ Cerec
- ♦ Primescan
- ♦ Medit

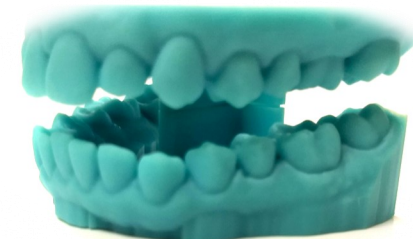

[dreamsystemsdigital@gmail.com](mailto:dreamsystemsdigital@gmail.com)



**myerson ema**  
...Sleep Better

510K—#K971794



### EMA ThermoFlex Hard

Made With Clasps  
Thermal Resin Compensates for  
Impression Defects & Missing Teeth  
Flat Plane Opposing Splint Design

### EMA VacuForm Splints

Standard with 2mm Splints  
Available in Hard / Soft Splints  
Requires Very Accurate  
PVS Impressions or Digital Models



### EMA First Step

90 Day Trial Device  
0.75 MM Vacuform Splints  
Three Yellow Strap Sizes 21,19,17mm  
Requires Very Accurate



### Advancement Strap Pack

Straps provided for 5mm of  
Advancement from Initial Setting  
5 Pairs Blue, 2 Pairs Yellow  
Bites Sent Advance of End to End  
will Receive 1-2 Pairs of Straps to  
Reduce Advancement, if Needed.

### Available Options

**Straps**  
**Sizes**  
21mm  
20mm  
19mm  
18mm  
17mm  
16mm  
15mm  
14mm  
13mm

White—Very Flexible  
(Provided Only Upon  
Request)  
Yellow—Flexible  
Blue—Standard  
Strap Less Flexible  
Clear—Rigid  
Provided Only Upon  
Request for Less  
Movement Side to Side



### Invisalign Retainers



Upper Invisalign Tray Used  
in place of Oasys Cushion




The Lower OASYS  
Splint is Designed to  
Function with Lower  
Tray Transitions



## Unique Solutions . . . for Challenging Dentition

### Minimal Dentition

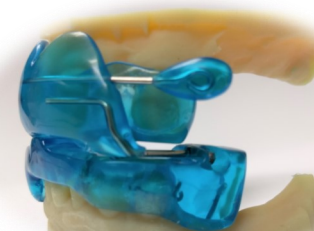

### Maxillary Denture





## Premium & Medicare OASYS Herbst Designs

### Premium Welded

Fully-Welded Frame, OASYS Herbst Hinges with Clasps, Ball Hooks & Elastics


### Medicare

OASYS Herbst Hinges with Welded Pivots; Clasps, Ball Hooks & Elastics



## 10 MM Advancement

MM Scale Shows Advancement

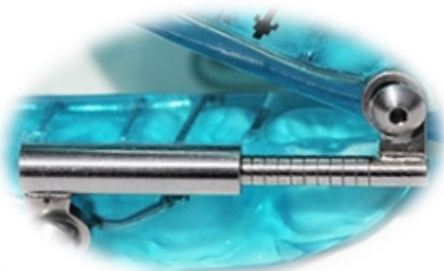

## Intra-Oral Activation

Easy Advancement

One Revolution = 1/2 MM

Two Revolutions = 1 MM

## Available Options

Medicare Compliant Nasal Dilators and Tongue Lifters  
Comfort Caps Available for Side Screw Irritation



**510K—#K083209**  
**Medicare Compliant**  
**#M001 OASYS Hinge Appliance**  
**+Nasal Dilators—#M002**  
**+Tongue Lifters—#M003**



## Premium & Medicare Herbst Designs

### Premium Welded

Fully-Welded Frame, Standard-Herbst Hinges, with Clasps, Ball Hooks & Elastics

### Medicare

Herbst Hinges with Welded Pivots; Clasps, Ball Hooks & Elastics



## 5 MM Advancement

Directional Arrows

Right Screw Up—Left Screw Down



## Additional Advancement

2 MM Crimpable Shims Available  
Longer Arms for Additional 5 MM  
Comfort Caps Available for Side Screw Irritation



**510K—#K083209**  
**Medicare Compliant—#MC1**  
**+Nasal Dilators—#M002**  
**+Tongue Lifters—#M003**

